



# Difa<sup>®</sup> Iodine

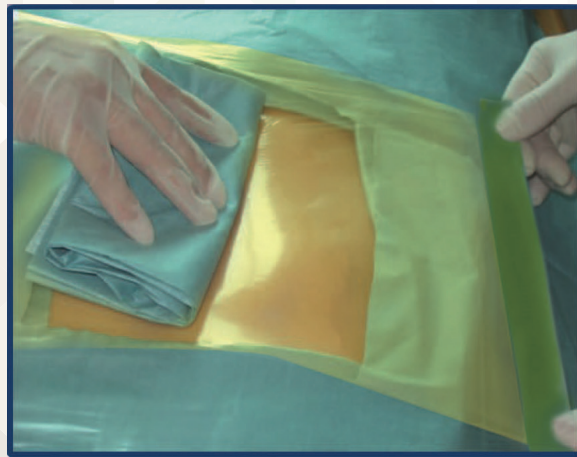

[www.zaditibbigerecler.com](http://www.zaditibbigerecler.com)

# Difa<sup>®</sup> Iodine

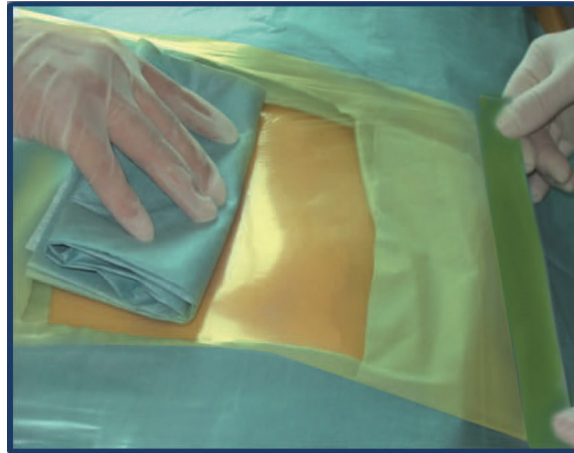

## Description

**Difa<sup>®</sup> Iodine –Theatre Incise Drape.**

Difa Iodine Incise Drape is a yellowish coloured, transparent, adhesive Polyurethane film which allows the skin to breathe, preventing moisture build-up under the drape. As a result adheres to the surrounding skin throughout long surgical operations and more importantly, right up to the wound edge. The adhesive used contains PVP-iodine, providing continuous broad spectrum antimicrobial activity, and a superior adhesion to the skin and wound edge.

**Difa<sup>®</sup> Iodine Frame Style** is the new generation of incise drapes which have a unique frame delivery system, enabling easy application by only one person, rather than two persons application in ordinary incise drapes.

## Features

- Frame border gives support to the edges of Difa Iodine Frame style, preventing wrinkles during application.
- Adheres right up to the wound edge.
- Yellowish Transparent surface.
- Thin but extremely strong.
- Conformable and elastic.
- Provides a sterile working field around the operation site and a broad spectrum antimicrobial activity.
- Holds the drapes securely in place.
- Sterile and individually wrapped in peel apart pouch.
- Sterilized by gamma irradiation.
- Bacterial proof and waterproof film.
- Available in a variety of sizes to meet your procedural and incision-size needs.

## Benefits

- **Adheres right up to wound edge**  
The drape provides protection from lateral infection because it stays in place even throughout major operations involving prolonged retraction of wound edges.
- **Transparent non-glare surface**  
The incision site is clearly visible without glare.
- **Thin but extremely strong**  
Difa Iodine Drape is tear resistant, but can be cleanly cut with a scalpel.
- **Conformable and elastic**  
Difa Iodine Drape is very conformable on operation site, its elasticity makes it suitable for use where heavy manipulative surgery is present, for example, in total hip replacement.
- **Maintains a sterile operation site**  
Difa Iodine Incise Drape provides a sterile working field, on which internal organs may be placed without fear of contamination.
- **Continuous broad spectrum antimicrobial activity**  
Iodine containing adhesive is a broad spectrum antimicrobial.

# Difa<sup>®</sup> Iodine

## Indications

**Difa iodine incise drapes:** are appropriate for surgical procedures where there is concern about bacterial wound contamination from skin flora and specifically surgeries such as:


- Orthopaedics.
- Neurosurgery.
- Major abdominal and implant surgery.
- Plastic Surgery.
- Open Heart and Thoracic.
- Ophthalmic.
- Paediatric.

## Contraindications:


Do not use Difa iodine on patient with known sensitivity to iodine

## Application:

### Difa Iodine

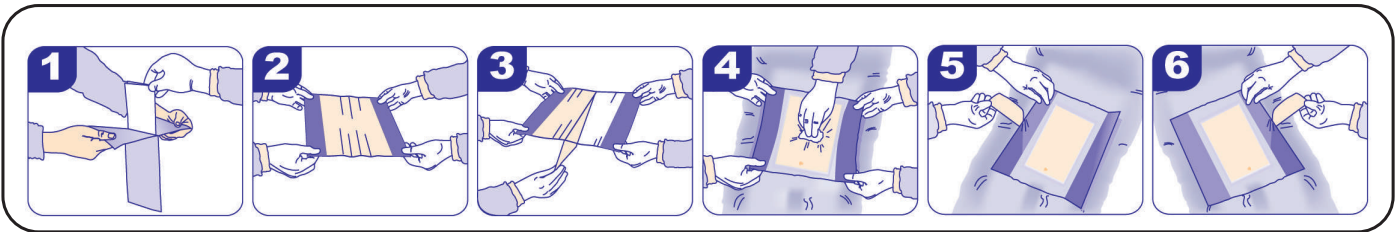

### Difa Iodine Frame Style

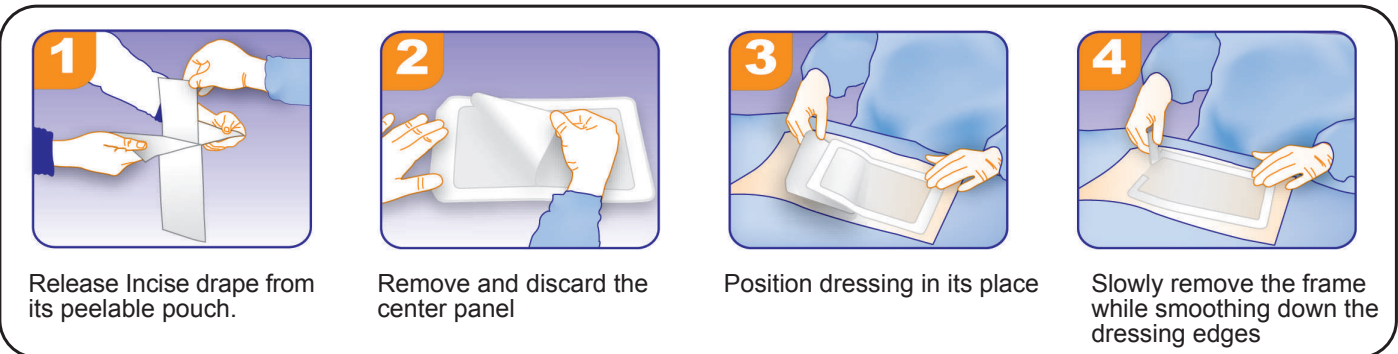

Product code	Size	Drapes per Case
Difa 1520 Iodine	15x20 cm	40
Difa 1528 Iodine	15x28 cm	40
Difa 3028 Iodine	30x28 cm	40
Difa 3030 Iodine	30x30 cm	40
Difa 2845 Iodine	28x45 cm	40
Difa 3045 Iodine	30x45cm	40
Difa 4555 Iodine	45x55 cm	40
Difa 4045 Iodine	40x45 cm	40
Difa 4245 Iodine	42x45 cm	40
Difa 3060 Iodine	30x60 cm	40
Difa 4560 Iodine	45x60 cm	40

Product code	Size	Drapes per Case
Difa 5080 Iodine	50x80 cm	40
Difa 9080 Iodine	90x80 cm	40
Difa 9070 Iodine	90x70 cm	40
Difa 1014 IodineFS	10X14 cm	1200
Difa1010 Iodine FS	10x10 cm	1200
Difa1015 Iodine FS	10x15 cm	1200
Difa 1020 IodineFS	10x20 cm	1200
Difa1025 Iodine FS	10x25 cm	600
Difa1030 Iodine FS	10x30 cm	600
Difa 1520 IodineFS	15x20 cm	200
Difa 2030 IodineFS	20x30 cm	100

Other sizes available upon request

**Zadi Tıbbi Gereçler San. ve Tic. A.Ş**  
İOSB Mah. İmsan Sanayi Sit. B Blok No:13  
Küçükçekmece / İSTANBUL  
Tel: +90 212 495 73 03 Fax: +90 212 495 73 04  
Email: import@zaditibbigerecler.com

[www.zaditibbigerecler.com](http://www.zaditibbigerecler.com)