

Fibrosol® Family Dressings



Description:

Fibrosol® Product Line is a sterile dressing used to cover acute and chronic wounds that are excreting pus or other fluids;

Fibrosol® Non-Adhesive is made of 100% CMC fibers, that are carded to form a nonwoven, reinforced with yarns used to stitch bond the product for support.

Fibrosol® Extra Non-Adhesive is made of 100% CMC fibers, that are carded with yarn forming a double layered nonwoven stitched together to give high strength and more integrity during dressing removal.

Fibrosol® Foam Gentle is composed of soft absorbent polyurethane foam laminated to a **Fibrosol® Non adhesive** layer placed as an island pad under a waterproof, virus proof and bacteria proof polyurethane film coated with silicone adhesive that allow gentle dressing removal without pain. The silicone adhesive layer is present as a perforated layer covering the whole dressing.

Fibrosol® Foam Non-Adhesive is composed of non-adherent soft absorbent polyurethane foam laminated to a **Fibrosol® Non adhesive** layer under a waterproof, virus proof and bacteria proof polyurethane film. Additional tape or other means of fixation should be used.


Features:

- Provides excellent absorption and retention capabilities for moderate to highly exuding wounds.
- Minimization of “dead space” where fluid may collect and bacteria can grow: the formation of a soft translucent gel provides conformability to uneven wound surfaces, which in turn provides soft and conformable adding to the wound surface to form an intimate contact to the patient comfort.
- Creates a moist environment that aids wound healing.
- Acts as a bacterial barrier that prevents outer wound infection. As in case of “**Fibrosol® Foam dressings**” the outer film layer provides a waterproof viral and bacterial barrier which protects the wound from external contaminants.
- Bacteriostatic activity.
- Reduction in localized maceration and excoriation in case of “**Fibrosol® Foam dressing**”.
- Helps reduce wound pain while the dressing is in situ and upon removal.
- Hemostatic effect.
- Balanced absorption force due to gel blocking: The structure of dressing fibers provides sufficient strength and “allows fluid to rapidly permeate and fully expand the fibers. The result is a coherent gel that resists wicking within the fibers and prevents wicking between the fibers by gel blocking.”
- Maintenance of moisture balance: **Fibrosol® dressings** feature, routes for moisture permeability and fiber gellation, ensures that absorption and retention properties are achieved.
- Balanced inflammatory response: The trapping of liquid and harmful components inside the dressing helps to prevent them from reaching the surface of the wound.
- Silicone adhesive wound contact layer provides secure, skin-friendly, gentle adhesion in case of “**Fibrosol® Foam Gentle dressing**”.

Intended Purpose:

A- Under the supervision of a healthcare professional, the **Fibrosol® Non-Adhesive** and may be used for the management of:

- Partial thickness (second degree) burns;
- Diabetic foot ulcers, leg ulcers, (venous stasis ulcers, arterial ulcers and leg ulcers of mixed etiology) and pressure ulcers/sores (partial, full thickness and stage IHV);
- Wounds as an effective barrier to bacterial penetration of the dressing as this may help reduce infection;
- Surgical wounds left to heal by secondary intention such as dehisced post-operative surgical incisions, donor sites, dermatological;
- Surgical wounds that heal by primary intention such as dermatological and surgical incisions (e.g. orthopedic and vascular);
- Traumatic wounds;
- Wounds that are prone to bleeding, such as wounds that have been mechanically or surgically debrided and donor sites;
- Oncology wounds with exudate, such as fungoides-cutaneous tumors, fungating carcinoma, cutaneous metastasis, Kaposi's sarcoma, and angiosarcoma;
- Management of painful wounds.

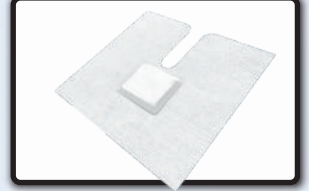

Fibrosol[®] Family Dressings

Sizes

Fibrosol [®] Non Adhesive	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIB5050	5 x 5 cm	10	48	
FIB100100	10 x 10 cm	10	48	
FIB150150	15 x 15 cm	10	24	
FIB20450	2 x 45 cm (ribbon)	5	24	
FIB25450	2.5 x 45 cm (ribbon)	5	24	
FIB100200	10 x 20 cm	10	24	
FIB100250	10 x 25 cm	5	24	
FIB100300	10 x 30 cm	5	24	
FIB150200	15 x 20 cm	5	24	
FIB150250	15 x 25 cm	5	24	
FIB20300	2 x 30 Cm	5	24	
FIB20400	2 x 40cm	5	24	
FIB200200	20 x 20 cm	5	24	
FIB200300	20 x 30 cm	5	24	
FIB50100	5 x 10 cm	10	48	
FIB50300	5 x 30 cm	5	24	
FIB5075	5 x 7.5 cm	10	48	
FIB70120	7 x 12 cm	10	48	
FIB75100	7.5 x 10 cm	10	48	
FIB75120	7.5 x 12 cm	10	48	
FIB75150	7.5 x 15cm	10	24	
FIB75200	7.5 x 20cm	10	24	
FIB5070	5 x 7cm	10	48	
FIB6070	6 x 7cm	10	48	
FIB100120	10 x 12cm	10	48	
FIB6080	6 x 8cm	10	48	
FIB165180	16.5 x 18cm	10	24	
FIB220220	22 x 22cm	10	24	
FIB7575	7.5 x 7.5 cm	10	48	
FIB125125	12,5x12,5	10	24	
FIB140150	14x15	10	24	
FIB175175	17,5x17,5	10	24	
FIB50400	5x40	10	24	
FIB100150	10x15 cm	10	24	
FIB135150	13,5x15	10	24	
FIB200240	20 x 24 cm	5	24	
FIB25400	2.5 x 40 cm (ribbon)	5	24	
Fibrosol [®] Extra Non Adhesive	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBEX5050	5 x 5 cm	10	48	
FIBEX100100	10 x 10 cm	10	48	
FIBEX150150	15 x 15 cm	10	24	
FIBEX20450	2 x 45 cm (ribbon)	5	24	
FIBEX25450	2.5 x 45 cm (ribbon)	5	24	
FIBEX100200	10 x 20 cm	10	24	
FIBEX100250	10 x 25 cm	10	24	
FIBEX100300	10 x 30 cm	10	24	
FIBEX150200	15 x 20 cm	10	24	
FIBEX150250	15 x 25 cm	5	24	

Fibrosol[®] Family Dressings



Fibrosol [®] Extra Non Adhesive	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBEX20300	2 x 30 Cm	5	24	
FIBEX20400	2 x 40cm	5	24	
FIBEX200200	20 x 20 cm	5	24	
FIBEX200300	20 x 30 cm	5	24	
FIBEX50100	5 x 10 cm	10	48	
FIBEX50300	5 x 30 cm	5	24	
FIBEX5075	5 x 7.5 cm	10	48	
FIBEX70120	7 x 12 cm	10	48	
FIBEX75100	7.5 x 10 cm	10	24	
FIBEX75120	7.5 x 12 cm	10	48	
FIBEX5070	5 x 7cm	10	48	
FIBEX6070	6 x 7cm	10	48	
FIBEX100120	10 x 12cm	10	48	
FIBEX6080	6 x 8cm	10	48	
FIBEX165180	16.5 x 18cm	10	24	
FIBEX220220	22 x 22cm	10	24	
FIBEX75150	7.5 x 15cm	10	24	
FIBEX75200	7.5 x 20cm	10	24	
FIBEX7575	7.5 x 7.5 cm	10	48	
FIBEX100150	10x15 cm	10	24	
FIBEX40100	4x10	10	48	
FIBEX40200	4x20	10	24	
FIBEX40300	4x30	10	24	
FIBEX50400	5x40	10	24	
FIBEX10450	1x45	5	24	
FIBEX125125	12,5x12,5	10	24	
FIBEX135150	13,5x15	10	24	
FIBEX180230	18x23	5	24	
FIBEX200240	20 x 24 cm	5	24	
FIBEX25400	2.5 x 40 cm (ribbon)	5	24	
Fibrosol [®] adhesive	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBAD5070	5x7cm	10	48	
FIBAD6070	6x7cm	10	48	
FIBAD6080	6x8cm	10	48	
FIBAD5090	5x9cm	10	48	
FIBAD6083	6x8.3cm	10	48	
FIBAD50100	5x10cm	10	48	
FIBAD6082	6x8.25cm	10	48	
FIBAD6085	6x8.5cm	10	48	
FIBAD7080	7x8cm	10	48	
FIBAD6580	6.5x8cm	10	48	
FIBAD7575	7.5x7.5cm	10	48	
FIBAD6090	6x9cm	10	48	
FIBAD60100	6x10cm	10	48	
FIBAD5082	5x8.25cm	10	48	
FIBAD75100	7.5x10cm	10	48	
FIBAD80100	8x10cm	10	48	
FIBAD82100	8.25x10cm	10	48	
FIBAD100100	10x10cm	10	48	
FIBAD70120	7x12cm	10	48	
FIBAD80120	8x12cm	10	48	

Fibrosol® Family Dressings

B- **Fibrosol® Foam Dressings; Fibrosol® Foam Gentle and Fibrosol® Foam Non-Adhesive** may be used for the management of both chronic and acute wounds, such as:

- Leg ulcers, pressure ulcers (Stage II-IV) and diabetic ulcers;
- Surgical wounds (post-operative, donor sites, dermatological);
- Partial thickness (second degree) burns;
- Traumatic or surgical wounds left to heal by secondary intention such as dehisced surgical incisions;
- Surgical wounds that heal by primary intent, such as dermatological and surgical incisions (e.g. orthopedic and vascular);
- Local management of wounds that are prone to bleeding, such as wounds that have been mechanically or surgically debrided and donor sites;
- Abrasions;
- Lacerations;
- Minor cuts;
- Painful wounds;
- Minor scalds and burns;
- **Fibrosol® Foam Dressings** may be included in a comprehensive protocol of care to protect against skin breakdown.

Then add table of sizes and codes for **Fibrosol® Non-Adhesive, Fibrosol® Extra Non-Adhesive, Fibrosol® Foam Gentle** and **Fibrosol® Foam Non-Adhesive Only**

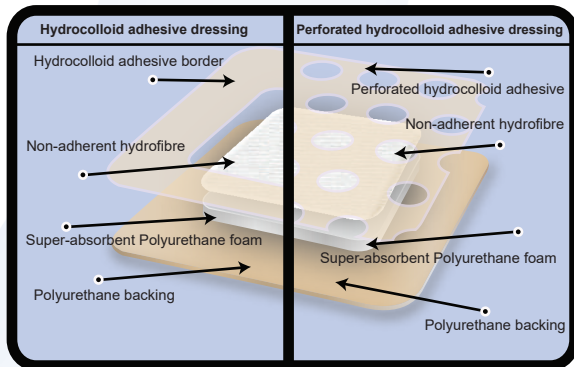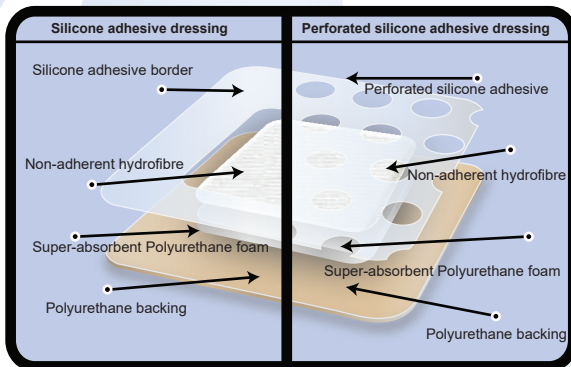

Fibrosol® Family Dressings

Sizes

Fibrosol® Non Adhesive	Fibrosol® Extra Non Adhesive	cm	inch	Pieces/Box	Boxes/Case
FIB5050	FIBEX5050	5 cm x 5 cm	2" x 2"	10	48
FIB1010	FIBEX1010	10 cm x 10 cm	4" x 4"	10	24
FIB1515	FIBEX1515	15 cm x 15 cm	6" x 6"	10	24
FIB245	FIBEX245	2 cm x 45 cm	3/4" x 18"	5	24
Fibrosol® adhesive	Fibrosol® Extra adhesive	cm	inch	Pieces/Box	Boxes/Case
FIBAD7050	FIBADEX7050	7cmx5cm	2¾ x 2 in	50	32
FIBAD1010	FIBADEX1010	10cmx10cm	4 x 4 in	40	32
FIBAD1015	FIBADEX1015	10cmx15cm	4 x 6 in	20	32
FIBAD1020	FIBADEX1020	10cmx20cm	4 x 8 in	20	32
FIBAD1025	FIBADEX1025	10cmx25cm	4 x 10 in	20	32
FIBAD1030	FIBADEX1030	10cmx30cm	4 x 12 in	20	32
Fibrosol® adhesive Film	Fibrosol® Extra adhesive Film	cm	inch	Pieces/Box	Boxes/Case
FIBPUAD7075	FIBPUADEX7050	7cmx5cm	2¾ x 2 in	50	32
FIBPUAD1010	FIBPUADEX1010	10cmx10cm	4 x 4 in	40	32
FIBPUAD1015	FIBPUADEX1015	10cmx15cm	4 x 6 in	20	32
FIBPUAD1020	FIBPUADEX1020	10cmx20cm	4 x 8 in	20	32
FIBPUAD1025	FIBPUADEX1025	10cmx25cm	4 x 10 in	20	32
FIBPUAD1030	FIBPUADEX1030	10cmx30cm	4 x 12 in	20	32
Fibrosol® adhesive Film Frame Style	Fibrosol® Extra adhesive Film Frame style	cm	inch	Pieces/Box	Boxes/Case
FIBPUFSAD7060	FIBPUFSADEX7060	7 cm x 6 cm	2¾ x 2¼ in	50	32
FIBPUFSAD1010	FIBPUFSADEX1010	10cmx10cm	4 x 4 in	40	32
FIBPUFSAD1015	FIBPUFSADEX1015	10cmx15cm	4 x 6 in	20	32
FIBPUFSAD1020	FIBPUFSADEX1020	10cmx20cm	4 x 8 in	20	32
FIBPUFSAD1025	FIBPUFSADEX1025	10cmx25cm	4 x 10 in	20	32
FIBPUFSAD1030	FIBPUFSADEX1030	10cmx30cm	4 x 12 in	20	32
Fibrosol® adhesive Ultra	Fibrosol® Extra adhesive Ultra	cm	inch	Pieces/Box	Boxes/Case
FIBULTRA7050	FIBULTRAEX7050	7cmx5cm	2¾ x 2 in	50	32
FIBULTRA1010	FIBULTRAEX1010	10cmx10cm	4 x 4 in	40	32
FIBULTRA1015	FIBULTRAEX1015	10cmx15cm	4 x 6 in	20	32
FIBULTRA1020	FIBULTRAEX1020	10cmx20cm	4 x 8 in	20	32
FIBULTRA1025	FIBULTRAEX1025	10cmx25cm	4 x 10 in	20	32
FIBULTRA1030	FIBULTRAEX1030	10cmx30cm	4 x 12 in	20	32
Fibrosol® IV	Fibrosol® Extra IV	cm	inch	Pieces/Box	Boxes/Case
FIBIV76	FIBIVEX76	7 cm x 6 cm	2¾ x 2¼ in	100	16
FIBIV68	FIBIVEX68	6 cm x 8 cm	4 x 4¾ in	100	16
FIBIV1012	FIBIVEX1012	10 cm x 12 cm	2¾ x 3½ in	100	16
Fibrosol® IV Film	Fibrosol® Extra IV Film	cm	inch	Pieces/Box	Boxes/Case
FIBIVPU76	FIBIVPUEX76	7 cm x 6 cm	2¾ x 2¼ in	100	16
FIBIVPU68	FIBIVPUEX68	6 cm x 8 cm	4 x 4¾ in	100	16
FIBIVPU1012	FIBIVPUEX1012	10 cm x 12 cm	2¾ x 3½ in	100	16
Fibrosol® Foam Gentle	Fibrosol® Extra Foam Gentle	cm	inch	Pieces/Box	Boxes/Case
FIBFOAMG5050	FIBFOAMGEX5050	5 cmx5 cm	2 x 2 in	10	48
FIBFOAMG1010	FIBFOAMGEX1010	10cmx10cm	4 x 4 in	5	24
FIBFOAMG1020	FIBFOAMGEX1020	10cmx20cm	4 x 8 in	5	24
FIBFOAMG1515	FIBFOAMGEX1515	15cmx15cm	6 x 6 in	5	24
FIBFOAMG2020	FIBFOAMGEX2020	20cmx20cm	8 x 8 in	5	24
FIBFOAMG16518	FIBFOAMGEX16518	18cmx16.5cm	7¼ x 6½	5	24
FIBFOAMG2222	FIBFOAMGEX2222	22cmx22cm	8 ⅝ x 8 ⅝	3	24

Fibrosol® Family Dressings

Fibrosol® Foam Basic Frame	Fibrosol® Extra Foam Basic Frame	cm	inch	Pieces/Box	Boxes/Case
FIBFOAMBF1010	FIBFOAMBFEX1010	10cmx10cm	4 x 4 in	10	48
FIBFOAMBF1020	FIBFOAMBFEX1020	10cmx20cm	4 x 8 in	5	24
FIBFOAMBF1515	FIBFOAMBFEX1515	15cmx15cm	6 x 6 in	5	24
FIBFOAMBF2020	FIBFOAMBFEX2020	20cmx20cm	8 x 8 in	3	24
FIBFOAMBF18165	FIBFOAMBFEX18165	18cmx16.5cm	7 ³ / ₄ x 6 ¹ / ₂	3	24
FIBFOAMBF2222	FIBFOAMBFEX2222	22cmx22cm	8 ⁵ / ₈ x 8 ⁵ / ₈	3	24
Fibrosol® Foam Basic Perforated	Fibrosol® Extra Foam Basic Perforated	cm	inch	Pieces/Box	Boxes/Case
FIBFOAMBP1010	FIBFOAMBPPEX1010	10cmx10cm	4 x 4 in	10	48
FIBFOAMBP1020	FIBFOAMBPPEX1020	10cmx20cm	4 x 8 in	5	24
FIBFOAMBP1515	FIBFOAMBPPEX1515	15cmx15cm	6 x 6 in	5	24
FIBFOAMBP2020	FIBFOAMBPPEX2020	20cmx20cm	8 x 8 in	3	24
FIBFOAMBP18165	FIBFOAMBPPEX18165	18cmx16.5cm	7 ³ / ₄ x 6 ¹ / ₂	3	24
FIBFOAMBP2222	FIBFOAMBPPEX2222	22cmx22cm	8 ⁵ / ₈ x 8 ⁵ / ₈	3	24
Fibrosol® OTC Bandage		cm	inch	Pieces/Box	Boxes/Case
FIBOTC938		0.9cmx3.8cm	3/ ₈ x 1 ¹ / ₂ in	200	60
FIBOTC1938		1.9cmx3.8 cm	3/ ₄ x 1 ¹ / ₂ in	100	60
FIBOTC1963		1.9cmx6.3cm	3/ ₄ x 2 ¹ / ₂ in	100	60
FIBOTC1972		1.9cmx7.2cm	3/ ₄ x 2 ⁷ / ₈ in	100	60
FIBOTC2572		2.5cmx7.2cm	1 x 2 ⁷ / ₈ in	100	60
FIBOTC25		spot 2.5 cm	1 in	100	192
FIBOTC2525		2.5cmx2.5cm	1 x 1 in	200	60
FIBOTC5072		5cmx7.2cm	2 x 2 ⁷ / ₈ in	50	60
FIBOTC5060 <small>Fingertip</small>		5 cm x 6 cm	2 x 2 ³ / ₈ in	50	60
FIBOTC4563 <small>Butterfly</small>		4.5cmx6.3cm	1 ³ / ₄ x 2 ¹ / ₂ in	50	60
FIBOTC510		5 cm x 10 cm	2 x 4 in	50	60
Fibrosol® Dressing strip		cm	inch	Pieces/Box	Boxes/Case
FIBDERMI610		6cmx10cm	2 ³ / ₈ x 4 in	10	120
FIBDERMI810		8 cmx10 cm	3 ¹ / ₈ x 4 in	10	120
FIBDERMI6100		6cmx100cm	2 ³ / ₈ x 40 in	1	120
FIBDERMI8100		8 cmx100 cm	3 ¹ / ₈ x 40 in	1	120

Other sizes available upon request

Fibrosol Adhesive dressings are available with:

- Acrylic Adhesive
- Hot melt Adhesive
- HydroColloid Adhesive
- Silicone Adhesive
- Frame Adhesive layer
- Perforated Adhesive layer
- Loose non-woven
- Stich non-woven

Fibrosol® Family Dressings

Fibrosol® adhesive	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBAD100120	10x12cm	10	48	
FIBAD83120	8.3x12cm	10	48	
FIBAD90150	9x15cm	10	48	
FIBAD80150	8x15cm	10	48	
FIBAD82150	8.25x15cm	10	48	
FIBAD100140	10x14cm	10	48	
FIBAD100150	10x15cm	10	48	
FIBAD150150	15x15cm	10	24	
FIBAD83180	8.3x18cm	10	48	
FIBAD82200	8.25x20cm	10	24	
FIBAD90200	9x20cm	10	24	
FIBAD100200	10x20cm	10	24	
FIBAD200200	20x20cm	10	24	
FIBAD83240	8.3x24cm	10	24	
FIBAD82250	8.25x25cm	10	24	
FIBAD90250	9x25cm	10	24	
FIBAD100250	10x25cm	10	24	
FIBAD82300	8.25x30cm	10	24	
FIBAD90300	9x30cm	10	24	
FIBAD100300	10x30cm	10	24	
FIBAD90350	9x35cm	10	24	
FIBAD82350	8.25x35cm	10	24	
FIBAD100350	10x35cm	10	24	
FIBAD83415	8.3x41.5cm	10	24	
FIBAD100160	10 x 16 cm	10	48	
FIBAD100210	10 x 21 cm	10	24	
FIBAD100260	10 x 26 cm	10	24	
FIBAD100400	10 x 40 cm	10	24	
FIBAD120160	12 x 16 cm	10	48	
FIBAD2572	2.5 x 7.2 cm	400	8	
FIBAD5072	5 x 7.2 cm	10	48	
FIBAD5075	5 x 7.5 cm	10	48	
FIBAD6075	6 x 7.5 cm	10	48	
FIBAD70100	7 x 10 cm	10	48	
FIBAD85155	8.5 x 15.5 cm	10	48	
FIBAD85150	8.5 x 15cm	10	48	
FIBAD8595	8.5 x 9.5 cm	10	48	
FIBAD90100	9 x 10cm	10	48	
Fibrosol® Extra adhesive	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBADEX5070	5x7cm	10	48	
FIBADEX6070	6x7cm	10	48	
FIBADEX6080	6x8cm	10	48	
FIBADEX5090	5x9cm	10	48	
FIBADEX6083	6x8.3cm	10	48	
FIBADEX50100	5x10cm	10	48	
FIBADEX6082	6x8.25cm	10	48	
FIBADEX6085	6x8.5cm	10	48	
FIBADEX7080	7x8cm	10	48	
FIBADEX6580	6.5x8cm	10	48	
FIBADEX7575	7.5x7.5cm	10	48	
FIBADEX6090	6x9cm	10	48	

Fibrosol[®] Family Dressings

Fibrosol [®] Extra adhesive	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBADEX60100	6x10cm	10	48	
FIBADEX5082	5x8.25cm	10	48	
FIBADEX75100	7.5x10cm	10	48	
FIBADEX80100	8x10cm	10	48	
		100	8	
FIBADEX82100	8.25x10cm	10	48	
FIBADEX100100	10x10cm	10	48	
FIBADEX70120	7x12cm	10	48	
FIBADEX80120	8x12cm	10	48	
FIBADEX100120	10x12cm	10	48	
FIBADEX83120	8.3x12cm	10	48	
FIBADEX90150	9x15cm	10	48	
FIBADEX80150	8x15cm	10	48	
FIBADEX82150	8.25x15cm	10	48	
FIBADEX100140	10x14cm	10	48	
FIBADEX100150	10x15cm	10	48	
FIBADEX150150	15x15cm	10	24	
FIBADEX83180	8.3x18cm	10	48	
FIBADEX82200	8.25x20cm	10	24	
FIBADEX90200	9x20cm	10	24	
FIBADEX100200	10x20cm	10	24	
FIBADEX200200	20x20cm	10	24	
FIBADEX83240	8.3x24cm	10	24	
FIBADEX82250	8.25x25cm	10	24	
FIBADEX90250	9x25cm	10	24	
FIBADEX100250	10x25cm	10	24	
FIBADEX82300	8.25x30cm	10	24	
FIBADEX90300	9x30cm	10	24	
FIBADEX100300	10x30cm	10	24	
FIBADEX90350	9x35cm	10	24	
FIBADEX82350	8.25x35cm	10	24	
FIBADEX100350	10x35cm	10	24	
FIBADEX83415	8.3x41.5cm	10	24	
FIBADEX100160	10 x 16 cm	10	48	
FIBADEX100210	10 x 21 cm	10	24	
FIBADEX100260	10 x 26 cm	10	24	
FIBADEX100400	10 x 40 cm	10	24	
FIBADEX120160	12 x 16 cm	10	48	
FIBADEX2572	2.5 x 7.2 cm	400	8	
FIBADEX5072	5 x 7.2 cm	10	48	
FIBADEX5075	5 x 7.5 cm	10	48	
FIBADEX6075	6x7.5cm	10	48	
		200	8	
FIBADEX70100	7 x 10 cm	10	48	
FIBADEX85155	8.5 x 15.5 cm	10	48	
FIBADEX85150	8.5 x 15cm	10	48	
FIBADEX8595	8.5 x 9.5 cm	10	48	
FIBADEX90100	9 x 10cm	10	48	
Fibrosol [®] adhesive Film	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBPUAD5070	5x7cm	10	48	
FIBPUAD6070	6x7cm	10	48	

Fibrosol[®] Family Dressings

Fibrosol [®] adhesive Film	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBPUAD6080	6x8cm	10	48	
FIBPUAD5090	5x9cm	10	48	
FIBPUAD6083	6x8.3cm	10	48	
FIBPUAD50100	5x10cm	10	48	
FIBPUAD6082	6x8.25cm	10	48	
FIBPUAD6085	6x8.5cm	10	48	
FIBPUAD7080	7x8cm	10	48	
FIBPUAD6580	6.5x8cm	10	48	
FIBPUAD7575	7.5x7.5cm	10	48	
FIBPUAD6090	6x9cm	10	48	
FIBPUAD60100	6x10cm	10	48	
FIBPUAD5082	5x8.25cm	10	48	
FIBPUAD75100	7.5x10cm	10	48	
FIBPUAD80100	8x10cm	10	48	
FIBPUAD82100	8.25x10cm	10	48	
FIBPUAD100100	10x10cm	10	48	
FIBPUAD70120	7x12cm	10	48	
FIBPUAD80120	8x12cm	10	48	
FIBPUAD100120	10x12cm	10	48	
FIBPUAD83120	8.3x12cm	10	48	
FIBPUAD90150	9x15cm	10	48	
FIBPUAD80150	8x15cm	10	48	
FIBPUAD82150	8.25x15cm	10	48	
FIBPUAD100140	10x14cm	10	48	
FIBPUAD100150	10x15cm	10	48	
FIBPUAD150150	15x15cm	10	24	
FIBPUAD83180	8.3x18cm	10	48	
FIBPUAD82200	8.25x20cm	10	24	
FIBPUAD90200	9x20cm	10	24	
FIBPUAD100200	10x20cm	10	24	
FIBPUAD200200	20x20cm	10	24	
FIBPUAD83240	8.3x24cm	10	24	
FIBPUAD82250	8.25x25cm	10	24	
FIBPUAD90250	9x25cm	10	24	
FIBPUAD100250	10x25cm	10	24	
FIBPUAD82300	8.25x30cm	10	24	
FIBPUAD90300	9x30cm	10	24	
FIBPUAD100300	10x30cm	10	24	
FIBPUAD90350	9x35cm	10	24	
FIBPUAD82350	8.25x35cm	10	24	
FIBPUAD100350	10x35cm	10	24	
FIBPUAD83415	8.3x41.5cm	10	24	
FIBPUAD100160	10 x 16 cm	10	48	
FIBPUAD100210	10 x 21 cm	10	24	
FIBPUAD100260	10 x 26 cm	10	24	
FIBPUAD100400	10 x 40 cm	10	24	
FIBPUAD120160	12 x 16 cm	10	48	
FIBPUAD2572	2.5 x 7.2 cm	400	8	
FIBPUAD5072	5 x 7.2 cm	10	48	
FIBPUAD5075	5 x 7.5 cm	10	48	
FIBPUAD6075	6 x 7.5 cm	10	48	
FIBPUAD70100	7 x 10 cm	10	48	

Fibrosol[®] Family Dressings

Fibrosol [®] adhesive Film	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBPUAD85155	8.5 x 15.5 cm	10	48	
FIBPUAD85150	8.5 x 15cm	10	48	
FIBPUAD8595	8.5 x 9.5 cm	10	48	
FIBPUAD90100	9 x 10cm	10	48	
Fibrosol [®] Extra adhesive Film	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBPUADEX6080	6x8cm	10	48	
FIBPUADEX5090	5x9cm	10	48	
FIBPUADEX6083	6x8.3cm	10	48	
FIBPUADEX50100	5x10cm	10	48	
FIBPUADEX6082	6x8.25cm	10	48	
FIBPUADEX6085	6x8.5cm	10	48	
FIBPUADEX7080	7x8cm	10	48	
FIBPUADEX6580	6.5x8cm	10	48	
FIBPUADEX7575	7.5x7.5cm	10	48	
FIBPUADEX6090	6x9cm	10	48	
FIBPUADEX60100	6x10cm	10	48	
FIBPUADEX5082	5x8.25cm	10	48	
FIBPUADEX75100	7.5x10cm	10	48	
FIBPUADEX80100	8x10cm	10	48	
FIBPUADEX82100	8.25x10cm	10	48	
FIBPUADEX100100	10x10cm	10	48	
FIBPUADEX70120	7x12cm	10	48	
FIBPUADEX80120	8x12cm	10	48	
FIBPUADEX100120	10x12cm	10	48	
FIBPUADEX83120	8.3x12cm	10	48	
FIBPUADEX90150	9x15cm	10	48	
FIBPUADEX80150	8x15cm	10	48	
FIBPUADEX82150	8.25x15cm	10	48	
FIBPUADEX100140	10x14cm	10	48	
FIBPUADEX100150	10x15cm	10	48	
FIBPUADEX150150	15x15cm	10	24	
FIBPUADEX83180	8.3x18cm	10	48	
FIBPUADEX82200	8.25x20cm	10	24	
FIBPUADEX90200	9x20cm	10	24	
FIBPUADEX100200	10x20cm	10	24	
FIBPUADEX200200	20x20cm	10	24	
FIBPUADEX83240	8.3x24cm	10	24	
FIBPUADEX82250	8.25x25cm	10	24	
FIBPUADEX90250	9x25cm	10	24	
FIBPUADEX100250	10x25cm	10	24	
FIBPUADEX82300	8.25x30cm	10	24	
FIBPUADEX90300	9x30cm	10	24	
FIBPUADEX100300	10x30cm	10	24	
FIBPUADEX90350	9x35cm	10	24	
FIBPUADEX82350	8.25x35cm	10	24	
FIBPUADEX100350	10x35cm	10	24	
FIBPUADEX83415	8.3x41.5cm	10	24	
FIBPUADEX100160	10 x 16 cm	10	48	
FIBPUADEX100210	10 x 21 cm	10	24	
FIBPUADEX100260	10 x 26 cm	10	24	
FIBPUADEX100400	10 x 40 cm	10	24	
FIBPUADEX120160	12 x 16 cm	10	48	

Fibrosol® Family Dressings

Fibrosol® Extra adhesive Film	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBPUADEX2572	2.5 x 7.2 cm	400	8	
FIBPUADEX5072	5 x 7.2 cm	10	48	
FIBPUADEX5075	5 x 7.5 cm	10	48	
FIBPUADEX6075	6 x 7.5 cm	10	48	
FIBPUADEX70100	7 x 10 cm	10	48	
FIBPUADEX85155	8.5 x 15.5 cm	10	48	
FIBPUADEX85150	8.5 x 15cm	10	48	
FIBPUADEX8595	8.5 x 9.5 cm	10	48	
FIBPUADEX90100	9 x 10cm	10	48	
Fibrosol® adhesive Film Frame Style	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBPUFSAD6080	6x8cm	10	48	
FIBPUFSAD5090	5x9cm	10	48	
FIBPUFSAD6083	6x8.3cm	10	48	
FIBPUFSAD50100	5x10cm	10	48	
FIBPUFSAD6082	6x8.25cm	10	48	
FIBPUFSAD6085	6x8.5cm	10	48	
FIBPUFSAD7080	7x8cm	10	48	
FIBPUFSAD6580	6.5x8cm	10	48	
FIBPUFSAD7575	7.5x7.5cm	10	48	
FIBPUFSAD6090	6x9cm	10	48	
FIBPUFSAD60100	6x10cm	10	48	
FIBPUFSAD5082	5x8.25cm	10	48	
FIBPUFSAD75100	7.5x10cm	10	48	
FIBPUFSAD80100	8x10cm	10	48	
FIBPUFSAD82100	8.25x10cm	10	48	
FIBPUFSAD100100	10x10cm	10	48	
FIBPUFSAD70120	7x12cm	10	48	
FIBPUFSAD80120	8x12cm	10	48	
FIBPUFSAD100120	10x12cm	10	48	
FIBPUFSAD83120	8.3x12cm	10	48	
FIBPUFSAD90150	9x15cm	10	48	
FIBPUFSAD80150	8x15cm	10	48	
FIBPUFSAD82150	8.25x15cm	10	48	
FIBPUFSAD100140	10x14cm	10	48	
FIBPUFSAD100150	10x15cm	10	48	
FIBPUFSAD150150	15x15cm	10	24	
FIBPUFSAD83180	8.3x18cm	10	48	
FIBPUFSAD82200	8.25x20cm	10	24	
FIBPUFSAD90200	9x20cm	10	24	
FIBPUFSAD100200	10x20cm	10	24	
FIBPUFSAD200200	20x20cm	10	24	
FIBPUFSAD83240	8.3x24cm	10	24	
FIBPUFSAD82250	8.25x25cm	10	24	
FIBPUFSAD90250	9x25cm	10	24	
FIBPUFSAD100250	10x25cm	10	24	
FIBPUFSAD82300	8.25x30cm	10	24	
FIBPUFSAD90300	9x30cm	10	24	
FIBPUFSAD100300	10x30cm	10	24	
FIBPUFSAD90350	9x35cm	10	24	
FIBPUFSAD82350	8.25x35cm	10	24	
FIBPUFSAD100350	10x35cm	10	24	
FIBPUFSAD83415	8.3x41.5cm	10	24	

Fibrosol[®] Family Dressings

Fibrosol [®] adhesive Film Frame Style	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBPUFSAD100160	10 x 16 cm	10	48	
FIBPUFSAD100210	10 x 21 cm	10	24	
FIBPUFSAD100260	10 x 26 cm	10	24	
FIBPUFSAD100400	10 x 40 cm	10	24	
FIBPUFSAD120160	12 x 16 cm	10	48	
FIBPUFSAD2572	2.5 x 7.2 cm	400	8	
FIBPUFSAD5072	5 x 7.2 cm	10	48	
FIBPUFSAD5075	5 x 7.5 cm	10	48	
FIBPUFSAD6075	6 x 7.5 cm	10	48	
FIBPUFSAD70100	7 x 10 cm	10	48	
FIBPUFSAD85155	8.5 x 15.5 cm	10	48	
FIBPUFSAD85150	8.5 x 15cm	10	48	
FIBPUFSAD8595	8.5 x 9.5 cm	10	48	
FIBPUFSAD90100	9 x 10cm	10	48	
Fibrosol [®] Extra adhesive Film Frame style	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBPUFSADEX6080	6x8cm	10	48	
FIBPUFSADEX5090	5x9cm	10	48	
FIBPUFSADEX6083	6x8.3cm	10	48	
FIBPUFSADEX50100	5x10cm	10	48	
FIBPUFSADEX6082	6x8.25cm	10	48	
FIBPUFSADEX6085	6x8.5cm	10	48	
FIBPUFSADEX7080	7x8cm	10	48	
FIBPUFSADEX6580	6.5x8cm	10	48	
FIBPUFSADEX7575	7.5x7.5cm	10	48	
FIBPUFSADEX6090	6x9cm	10	48	
FIBPUFSADEX60100	6x10cm	10	48	
FIBPUFSADEX5082	5x8.25cm	10	48	
FIBPUFSADEX75100	7.5x10cm	10	48	
FIBPUFSADEX80100	8x10cm	10	48	
FIBPUFSADEX82100	8.25x10cm	10	48	
FIBPUFSADEX100100	10x10cm	10	48	
FIBPUFSADEX70120	7x12cm	10	48	
FIBPUFSADEX80120	8x12cm	10	48	
FIBPUFSADEX100120	10x12cm	10	48	
FIBPUFSADEX83120	8.3x12cm	10	48	
FIBPUFSADEX90150	9x15cm	10	48	
FIBPUFSADEX80150	8x15cm	10	48	
FIBPUFSADEX82150	8.25x15cm	10	48	
FIBPUFSADEX100140	10x14cm	10	48	
FIBPUFSADEX100150	10x15cm	10	48	
FIBPUFSADEX150150	15x15cm	10	24	
FIBPUFSADEX83180	8.3x18cm	10	48	
FIBPUFSADEX82200	8.25x20cm	10	24	
FIBPUFSADEX90200	9x20cm	10	24	
FIBPUFSADEX100200	10x20cm	10	24	
FIBPUFSADEX200200	20x20cm	10	24	
FIBPUFSADEX83240	8.3x24cm	10	24	
FIBPUFSADEX82250	8.25x25cm	10	24	
FIBPUFSADEX90250	9x25cm	10	24	
FIBPUFSADEX100250	10x25cm	10	24	
FIBPUFSADEX82300	8.25x30cm	10	24	
FIBPUFSADEX90300	9x30cm	10	24	

- Amria - Alexandria - Egypt
 Tel. : 002 03 4500264 Fax : 002 03 4500263
 - Kafr-El-Zayat P.O.Box 38 Egypt
 Tel. : 002 040 2545643 Fax : 002 040 2545645
 E-mail : sales@pharmaplast-online.com

Fibrosol[®] Family Dressings

Fibrosol [®] Extra adhesive Film Frame style	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBPUFSADEX100300	10x30cm	10	24	
FIBPUFSADEX90350	9x35cm	10	24	
FIBPUFSADEX82350	8.25x35cm	10	24	
FIBPUFSADEX100350	10x35cm	10	24	
FIBPUFSADEX83415	8.3x41.5cm	10	24	
FIBPUFSADEX100160	10 x 16 cm	10	48	
FIBPUFSADEX100210	10 x 21 cm	10	24	
FIBPUFSADEX100260	10 x 26 cm	10	24	
FIBPUFSADEX100400	10 x 40 cm	10	24	
FIBPUFSADEX120160	12 x 16 cm	10	48	
FIBPUFSADEX2572	2.5 x 7.2 cm	400	8	
FIBPUFSADEX5072	5 x 7.2 cm	10	48	
FIBPUFSADEX5075	5 x 7.5 cm	10	48	
FIBPUFSADEX6075	6 x 7.5 cm	10	48	
FIBPUFSADEX70100	7 x 10 cm	10	48	
FIBPUFSADEX85155	8.5 x 15.5 cm	10	48	
FIBPUFSADEX85150	8.5 x 15cm	10	48	
FIBPUFSADEX8595	8.5 x 9.5 cm	10	48	
FIBPUFSADEX90100	9 x 10cm	10	48	
Fibrosol [®] adhesive Ultra	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBULTRA5070	5x7cm	10	48	
FIBULTRA6070	6x7cm	10	48	
FIBULTRA6080	6x8cm	10	48	
FIBULTRA5090	5x9cm	10	48	
FIBULTRA6083	6x8.3cm	10	48	
FIBULTRA50100	5x10cm	10	48	
FIBULTRA6082	6x8.25cm	10	48	
FIBULTRA6085	6x8.5cm	10	48	
FIBULTRA7080	7x8cm	10	48	
FIBULTRA6580	6.5x8cm	10	48	
FIBULTRA7575	7.5x7.5cm	10	48	
FIBULTRA6090	6x9cm	10	48	
FIBULTRA60100	6x10cm	10	48	
FIBULTRA5082	5x8.25cm	10	48	
FIBULTRA75100	7.5x10cm	10	48	
FIBULTRA80100	8x10cm	10	48	
		100	8	
FIBULTRA82100	8.25x10cm	10	48	
FIBULTRA100100	10x10cm	10	48	
FIBULTRA70120	7x12cm	10	48	
FIBULTRA80120	8x12cm	10	48	
FIBULTRA100120	10x12cm	10	48	
FIBULTRA83120	8.3x12cm	10	48	
FIBULTRA90150	9x15cm	10	48	
FIBULTRA80150	8x15cm	10	48	
FIBULTRA82150	8.25x15cm	10	48	
FIBULTRA100140	10x14cm	10	48	
FIBULTRA100150	10x15cm	10	48	
FIBULTRA150150	15x15cm	10	24	
FIBULTRA83180	8.3x18cm	10	48	
FIBULTRA82200	8.25x20cm	10	24	
FIBULTRA90200	9x20cm	10	24	

Fibrosol[®] Family Dressings

Fibrosol [®] adhesive Ultra	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBULTRA100200	10x20cm	10	24	
FIBULTRA200200	20x20cm	10	24	
FIBULTRA83240	8.3x24cm	10	24	
FIBULTRA82250	8.25x25cm	10	24	
FIBULTRA90250	9x25cm	10	24	
FIBULTRA100250	10x25cm	10	24	
FIBULTRA82300	8.25x30cm	10	24	
FIBULTRA90300	9x30cm	10	24	
FIBULTRA100300	10x30cm	10	24	
FIBULTRA90350	9x35cm	10	24	
FIBULTRA82350	8.25x35cm	10	24	
FIBULTRA100350	10x35cm	10	24	
FIBULTRA83415	8.3x41.5cm	10	24	
FIBULTRA100160	10 x 16 cm	10	48	
FIBULTRA100210	10 x 21 cm	10	24	
FIBULTRA100260	10 x 26 cm	10	24	
FIBULTRA100400	10 x 40 cm	10	24	
FIBULTRA120160	12 x 16 cm	10	48	
FIBULTRA2572	2.5 x 7.2 cm	10	8	
FIBULTRA5072	5 x 7.2 cm	10	8	
FIBULTRA5075	5 x 7.5 cm	10	8	
FIBULTRA6075	6 x 7.5 cm	10	48	
		200	8	
FIBULTRA70100	7 x 10 cm	10	48	
FIBULTRA85155	8.5 x 15.5 cm	10	48	
FIBULTRA85150	8.5 x 15cm	10	48	
FIBULTRA8595	8.5 x 9.5 cm	10	48	
FIBULTRA90100	9 x 10cm	10	48	
Fibrosol [®] Extra adhesive Ultra	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBULTRAEX5070	5x7cm	10	48	
FIBULTRAEX60100	6x10cm	10	48	
FIBULTRAEX100100	10x10cm	10	48	
FIBULTRAEX100150	10x15cm	10	48	
FIBULTRAEX100200	10x20cm	10	24	
FIBULTRAEX100250	10x25cm	10	24	
FIBULTRAEX100300	10x30cm	10	24	
FIBULTRAEX100350	10x35cm	10	24	
FIBULTRAEX150150	15x15cm	10	24	
FIBULTRAEX200200	20x20cm	10	24	
FIBULTRAEX100160	10 x 16 cm	10	48	
FIBULTRAEX100210	10 x 21 cm	10	24	
FIBULTRAEX100260	10 x 26 cm	10	24	
FIBULTRAEX100400	10 x 40 cm	10	24	
FIBULTRAEX120160	12 x 16 cm	10	48	
FIBULTRAEX2572	2.5 x 7.2 cm	400	8	
FIBULTRAEX5072	5 x 7.2 cm	10	48	
FIBULTRAEX5075	5 x 7.5 cm	10	48	
FIBULTRAEX6075	6 x 7.5 cm	10	48	
FIBULTRAEX6085	6 x 8.5 cm	10	48	
FIBULTRAEX70100	7 x 10 cm	10	48	
FIBULTRAEX80100	8 x 10 cm	10	48	
FIBULTRAEX85155	8.5 x 15.5 cm	10	48	

Fibrosol[®] Family Dressings

Fibrosol [®] Extra adhesive Ultra	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBULTRAEX85150	8.5 x 15cm	10	48	
FIBULTRAEX8595	8.5 x 9.5 cm	10	48	
FIBULTRAEX90100	9 x 10cm	10	48	
Fibrosol [®] IV	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBIV6070	6X7CM	100	16	
FIBIV7080	7X8CM	100	16	
FIBIV82100	8.25X10CM	60	16	
FIBIV7085	7X8.5CM	60	16	
FIBIV7090	7X9CM	60	16	
FIBIV100120	10X12CM	60	16	
FIBIV6080	6X8CM	100	16	
FIBIV5057	5X5.7CM	100	16	
FIBIV7070	7x7 cm	100	16	
FIBIV80100	8x10 cm	100	16	
FIBIV6085	6x8.5 cm	100	16	
FIBIV85100	8.5x10 cm	100	16	
FIBIV85115	8.5x11.5 cm	100	16	
FIBIV90100	9x10 cm	100	16	
FIBIV100140	10x14 cm	60	16	
FIBIV100155	10x15.5 cm	60	16	
FIBIV5055	5x5.5 cm	100	16	
FIBIV6570	6.5x7 cm	100	16	
Fibrosol [®] Extra IV	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBIVEX7060	7X6CM	100	16	
FIBIVEX8070	8X7CM	100	16	
FIBIVEX82510	8.25X10CM	60	16	
FIBIVEX8570	8.5X7CM	60	16	
FIBIVEX97	9X7CM	60	16	
FIBIVEX1012	10X12CM	60	16	
FIBIVEX6080	6X8CM	100	16	
FIBIVEX5057	5X5.7CM	100	16	
FIBIVEX77	7x7 cm	100	16	
FIBIVEX810	8x10 cm	100	16	
FIBIVEX856	8.5x6 cm	100	16	
FIBIVEX8510	8.5x10 cm	100	16	
FIBIVEX85115	8.5x11.5 cm	100	16	
FIBIVEX910	9x10 cm	100	16	
FIBIVEX1014	10x14 cm	60	16	
FIBIVEX10155	10x15.5 cm	60	16	
FIBIVEX555	5.5x5 cm	100	16	
FIBIVEX657	6.5x7 cm	100	16	
Fibrosol [®] IV Film	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBIVPU6070	6X7CM	100	16	
FIBIVPU7070	7X7CM	100	16	
FIBIVPU7080	7X8CM	100	16	
FIBIVPU80100	8X10CM	60	16	
FIBIVPU82100	8.25X10CM	60	16	
FIBIVPU6085	6X8.5CM	100	16	
FIBIVPU7085	7X8.5CM	60	16	
FIBIVPU85100	8.5X10CM	60	16	
FIBIVPU7090	7X9CM	60	16	
FIBIVPU90100	9X10CM	60	16	

Fibrosol® Family Dressings

Fibrosol® IV Film	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBIVPU100120	10X12CM	60	16	
FIBIVPU100140	10X14CM	60	16	
FIBIVPU6080	6X8CM	100	16	
FIBIVPU85115	8.5X11.5CM	60	16	
FIBIVPU5057	5X5.7CM	100	16	
Fibrosol® Extra IV Film	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBIVPUEX6070	6X7CM	100	16	
FIBIVPUEX7070	7X7CM	100	16	
FIBIVPUEX7080	7X8CM	100	16	
FIBIVPUEX80100	8X10CM	60	16	
FIBIVPUEX82100	8.25X10CM	60	16	
FIBIVPUEX6085	6X8.5CM	100	16	
FIBIVPUEX7085	7X8.5CM	60	16	
FIBIVPUEX85100	8.5X10CM	60	16	
FIBIVPUEX7090	7X9CM	60	16	
FIBIVPUEX90100	9X10CM	60	16	
FIBIVPUEX100120	10X12CM	60	16	
FIBIVPUEX100140	10X14CM	60	16	
FIBIVPUEX6080	6X8CM	100	16	
FIBIVPUEX85115	8.5X11.5CM	60	16	
FIBIVPUEX5057	5X5.7CM	100	16	
Fibrosol® Foam Gentle	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBFOAMG75100	7.5x10cm	10	48	
FIBFOAMG100100	10x10cm	10	48	
FIBFOAMG100150	10x15cm	5	24	
FIBFOAMG100200	10x20cm	5	24	
FIBFOAMG100350	10x35 cm	5	24	
FIBFOAMG100250	10x25cm	5	24	
FIBFOAMG150150	15x15cm	5	24	
FIBFOAMG200200	20x20cm	3	24	
FIBFOAMG165180	Concave	3	24	
FIBFOAMG220220	Sacral	3	24	
FIBFOAMG125125	12.5x12.5cm	5	24	
FIBFOAMG100300	10x30cm	5	24	
FIBFOAMG170170	170x170 mm	3	24	
	Small Sacral			
FIBFOAMG5050	50x50 mm	10	48	
FIBFOAMG6060	60x60 mm	10	48	
FIBFOAMG100125	100x125mm	5	24	
FIBFOAMG175175	175x175 mm	3	24	
FIBFOAMG7575	7.5x7.5cm	10	24	
FIBFOAMG150200	15x20cm	3	24	
FIBFOAMG100120	10x12cm			
FIBFOAMG225225	22.5x22.5 cm			
FIBFOAMG150175	15x17.5			
FIBFOAMG175200	17.5x20			
FIBFOAMG250300	25x30			
FIBFOAMG5070	5x7			
FIBFOAMG170170N	170x170 mm			
Fibrosol® Extra Foam Gentle	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBFOAMGEX75100	7.5x10cm	10	48	
FIBFOAMGEX100100	10x10cm	10	48	

Fibrosol® Family Dressings

Fibrosol® Extra Foam Gentle	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBFOAMGEX100150	10x15cm	5	24	
FIBFOAMGEX100200	10x20 cm	5	24	
FIBFOAMGEX100350	10x35 cm	5	24	
FIBFOAMGEX100250	10x25cm	5	24	
FIBFOAMGEX150150	15x15cm	5	24	
FIBFOAMGEX200200	20x20cm	3	24	
FIBFOAMGEX165180	Concave	3	24	
FIBFOAMGEX220220	Sacral	3	24	
FIBFOAMGEX125125	12.5x12.5cm	5	24	
FIBFOAMGEX100300	10x30cm	5	24	
FIBFOAMGEX170170	170x170 mm	3	24	
	Small Sacral			
FIBFOAMGEX5050	50x50 mm	10	48	
FIBFOAMGEX6060	60x60 mm	10	48	
FIBFOAMGEX100125	100x125mm	5	24	
FIBFOAMGEX175175	175x175 mm	3	24	
FIBFOAMGEX7575	7.5x7.5cm	10	24	
FIBFOAMGEX150200	15x20cm	3	24	
FIBFOAMGEX170170N	170x170 mm			
Fibrosol® Foam Basic	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBFOAMB5050	50x50 mm	10	48	
FIBFOAMB7575	75x75 mm	10	48	
FIBFOAMB75100	75 x 100 mm	10	48	
FIBFOAMB100100	100 x 100 mm	10	48	
FIBFOAMB100150	100x150 mm	5	24	
FIBFOAMB150150	150 x 150 mm	5	24	
FIBFOAMB100200	100x200 mm	5	24	
FIBFOAMB100250	100x250 mm	5	24	
FIBFOAMB100300	100x300 mm	5	24	
FIBFOAMB150200	150x200 mm	3	24	
FIBFOAMB100350	100x350 mm	5	24	
FIBFOAMB200200	200 x 200 mm	3	24	
FIBFOAMB6060	60x60 mm	10	48	
FIBFOAMB100125	100x125 mm	5	24	
FIBFOAMB125125	125x125 mm	5	24	
FIBFOAMB170170	170x170 mm	3	24	
FIBFOAMB175175	175x175 mm	3	24	
FIBFOAMB165180	Concave 165 x 180 mm	3	24	
FIBFOAMB220220	Sacral 220 x 220 mm	3	24	
Fibrosol® Extra Foam Basic	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBFOAMBEX5050	50x50 mm	10	48	
FIBFOAMBEX7575	75x75 mm	10	48	
FIBFOAMBEX75100	75 x 100 mm	10	48	
FIBFOAMBEX100100	100 x 100 mm	10	48	
FIBFOAMBEX100150	100x150 mm	5	24	
FIBFOAMBEX150150	150 x 150 mm	5	24	
FIBFOAMBEX100200	100x200 mm	5	24	
FIBFOAMBEX100250	100x250 mm	5	24	
FIBFOAMBEX100300	100x300 mm	5	24	
FIBFOAMBEX150200	150x200 mm	3	24	
FIBFOAMBEX100350	100x350 mm	5	24	
FIBFOAMBEX200200	200 x 200 mm	3	24	

Fibrosol® Family Dressings

Fibrosol® Extra Foam Basic	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBFOAMBEX6060	60x60 mm	10	48	
FIBFOAMBEX100125	100x125 mm	5	24	
FIBFOAMBEX125125	125x125 mm	5	24	
FIBFOAMBEX170170	170x170 mm	3	24	
FIBFOAMBEX175175	175x175 mm	3	24	
FIBFOAMBEX165180	Concave 165 x 180 mm	3	24	
FIBFOAMBEX220220	Sacral 220 x 220 mm	3	24	
Fibrosol® OTC Bandage	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBOTC938	9 X 38 mm	100	60	
FIBOTC1657	16 X 57 mm	100	60	
FIBOTC1938	19 X 38 mm	100	60	
FIBOTC1963	19 X 63 mm	100	60	
FIBOTC1964	19 X 64 mm	100	60	
FIBOTC1972	19 X 72 mm	100	60	
FIBOTC1976	19 X 76 mm	100	60	
FIBOTC2572	25 X 72 mm	100	60	
FIBOTC2576	25 X 76 mm	100	60	
FIBOTC3055	30 X 55 mm	100	60	
FIBOTC3064	30 X 64 mm	100	60	
FIBOTC3072	30 X 72 mm	50	60	
FIBOTCANCHOR	38 X 72 mm Anchor			
FIBOTC3838	38 X 38 mm	100	60	
FIBOTC3872	38 X 72 mm	50	60	
FIBOTC5072	50 X 72 mm	50	60	
FIBOTC50100	50 X 100 mm	50	60	
FIBOTC5060	50 X 60 mm Fingertip	50	60	
FIBOTC4563	45 X 63 mm Butterfly	50	60	
FIBOTC25	SPOT 25 mm	100	192	
FIBOTC22	SPOT 22 mm	100	192	
FIBOTC2525	25 X 25 mm	100	192	
FIBOTCEYEPAD5062	Eyepad child	20	40	
FIBOTCEYEPAD5882	Eyepad Adult	20	40	
FIBOTCEYEPAD64110	ocular	20	40	
FIBOTCFEET	Feet	10	48	
FIBOTC2235	22 x 35 mm	100	60	
FIBOTC2665	26 x 65 mm	100	60	
FIBOTC3070	30 x 70 mm	100	60	
FIBOTC64110	64 X 110 mm	50	60	
FIBOTC5284	52 X 84 mm	50	60	
FIBOTC20120	20 x 120 mm	100	60	
FIBOTC30120	30 x 120 mm	100	60	
FIBOTC20180	20 x 180 mm	100	60	
FIBOTC30180	30 x 180 mm	100	60	
FIBOTC2582	25 x 82 mm	100	60	
FIBOTC2956	29 x 56 mm	100	60	
FIBOTCCHARACTERS	characters	100	60	
FIBOTC1955	19 x 55 mm	100	60	
FIBOTC4550	45 x 50 mm Butterfly	50	60	
Fibrosol® Dressing strip	cm	Pieces/Box	Boxes/Case	Case Dimensions
FIBDERMI60100	60 x 100 mm	10	120	
FIBDERMI601000	60 x 1000 mm	1	120	
FIBDERMI80100	80 X 100 mm	10	120	
FIBDERMI801000	80 x 1000 mm	1	120	

- Amria - Alexandria - Egypt
 Tel. : 002 03 4500264 Fax : 002 03 4500263
 - Kafr-El-Zayat P.O.Box 38 Egypt
 Tel. : 002 040 2545643 Fax : 002 040 2545645
 E-mail : sales@pharmaplast-online.com

 **Pharmaplast**
www.pharmaplast-online.com

Pharmaplast Amria 23512, Alexandria - Egypt www.pharmaplast-online.com		Brochure artwork template no: I-TM /AR/001-6 (Version 3)						
Customer : Pharmaplast		Product name : Fibrosol Family						
Issue date : 13-05-2015		Prepared by : Nader Nabil zaki						
Colors : CMYK		Dimension : A4						
Issue Number	Revision Date	Description of Change	Rationale For Change	Amended by	Effective date	Approved by		
						Artwork manager	Marketing department	RA Department
1	N/A	New Issue	N/A					
2	20/02/2019	change codes at table of size	to be identical with sap codes					
3	07/07/2019	- amend text	- Text aligned with CER					
4	02-02-2022	change codes at table of size	File Maintenance					
4	04-04-2022	text updated	File Maintenance					